





Thermal Resistance & Spectroradio meter System for LED (TRS-1000)

Brochure

Global Office of Lisun Electronics Inc.

http://www.Lisungroup.com

Lisun Group (Hong Kong) Limited

Add: Room 803, Chevalier House, 45-51 Chatham Road South, Tsim Sha Tsui, KL, HK

Tel: 00852-68852050 Fax: 00852-30785638

Email: SalesHK@Lisungroup.com

Lisun Electronics (Shanghai) Co., Ltd

Add: 113-114, No. 1 Building, Nanxiang Zhidi Industry Park, No. 1101, Huyi Road, Jiading

District, Shanghai, 201802, China

Tel: +86(21)5108 3341 Fax: +86(21)5108 3342

Email: SalesSH@Lisungroup.com

Lisun Electronics Inc. (USA)

Add: 445 S. Figueroa Street, Los Angeless, CA 90071, U.S.A.

Email: Sales@Lisungroup.com

Lisun China Factory

Add: NO. 37, Xiangyuan Road, Hangzhou City, Zhejiang Province, China

Tel: +86-189-1799-6096

Email: Engineering@Lisungroup.com

Leader in Lighting & Electrical Test Instruments

Rev. 10/21/2019



System configuration:

- A. LED Fast Scan Spectroradiometer
- B. T3 LED Thermal and Electrical Performance Analyzer
- C. Temperature Control Device
- D. Standard Light source
- E. 0.3m Integrating sphere

Brief Introduction:

Thermal parameter is the key index to evaluate the quality of LED. Different temperature will have significant different effects on <u>optical parameters</u>, <u>color parameters</u>, electricity parameters and life of LED. Therefore, design and control of the thermal dissipation is the most important thing for the manufacturers and users. Referring to the relative international standard, TRS1000 could conduct integration testing on optic, color, electricity and thermal; meanwhile record the changing curve of the parameters at real time.

Technical Specification:

A. Thermal parameters:

Junction temperature, thermal resistance and K factor

• Thermal sink range: 5°C ~ 90°C

Temperature control system accuracy: ±1℃

Junction temperature accuracy: ±1℃

- B. Optic and colorimetric Parameters: <u>Spectral power distribution</u>, chromaticity coordinate, correlative color temperature, <u>peak wavelength</u>, dominant wavelength, semi breadth, red ratio, rendering index, <u>luminous flux</u> and optical power
 - Wavelength range: 380~780nmWavelength accuracy: ±0.3nm
 - Repeatability of chromaticity coordinate: ±0.001 (Under standard illuminator A)
 - Sampling interval: 1nm/5nm
 Photometry repeatability: 1%

• CCT range: 1000~100000K

Dominant wavelength range: 380~700nm

Sensor: 2048 unit array CCDIntegration time: 6µs~10s

C. Electrical parameters: Forward voltage (Vf) and forward current(If)

Output current range: 1-3000mA

Output voltage range: 45V

Max LED power dissipation of LED: 20W

Clamp:

- A. If test thermal VS optical parameters and color parameters, the user will need a temperature controlled clamp to mounted the LED, the temperature controlled clamp can be designed according to users' products.
- B. If test optical parameters and color parameters, we have 4 clamps to mounted the LED.







Features:

- Record the changing curve of parameters at real time.
- Optical fiber connector: SMA 905
- Communication port: USB2.0
- Software already built in database function, enable to output Excel format file.
- Programmable set thermal sink temperature, and If of LED for the purpose of imitating different working situation of LED
- Test condition complies with corresponding international standard
- Operation system: Windows 2000/XP/VISTA

Working condition requirement:

Temperature range: 15°C ~35°C

Max. Humidity: 70% RH

Note:


Standard Light source:

OSRAM Standard Lamp to calibrate the spectrum and <u>luminous flux</u> with Lisun Lab certification. The data can be traced NIM. The Standard Lamp Source is used to calibrate the integrating sphere system. The different size of Integrating Sphere should choose the right power of standard lamp source


Integrating sphere

- Painting of integrating spheres is according to CIE Pub.No.84(1989)
- Material: Pure barium soleplate (BaSO4)
- BaSO4 coating: $\rho(\lambda) \ge 0.96(450 \text{nm} \sim 800 \text{nm})$ and $\rho(\lambda) \ge 0.92(380 \text{nm} \sim 450 \text{nm})$
- Fine diffuse reflection: Reflectancep \approx 0.8 and accuracy of p(λ)<1.5%


Software:

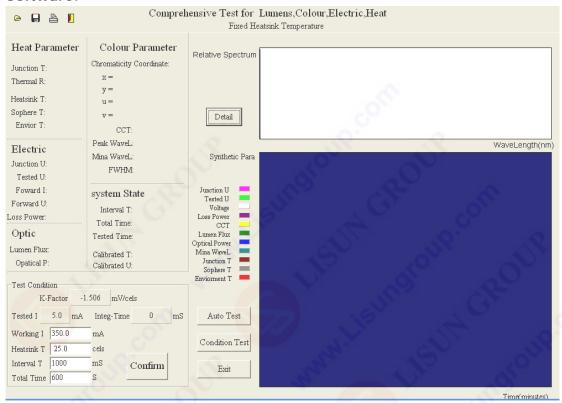


Comprehensive Test for Lumen, Color, Electric and Heat C

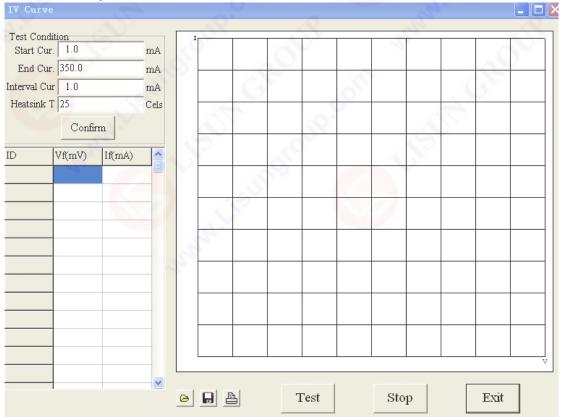

IV Curve Test

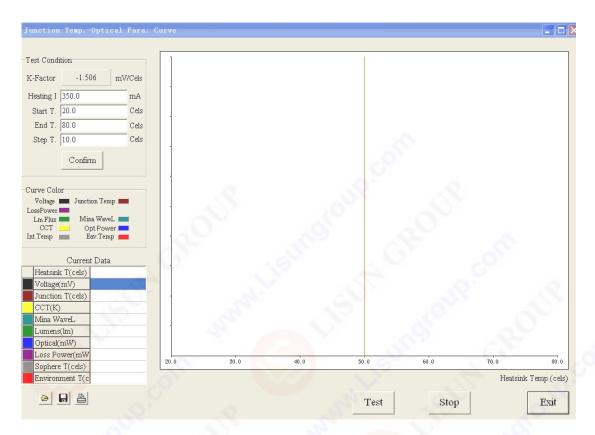

Junction Temp and Optical Curve

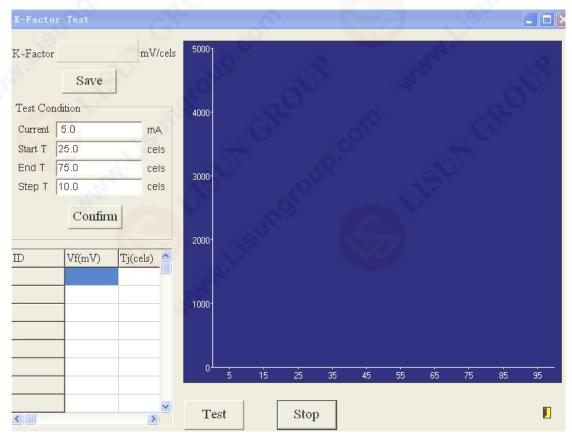

K Factor Test